

Company Profile

Starting in 1999, Blais Microscope Company initially only serviced the Midwest. We now service over 5000 accounts across 35 states, and we continue to grow. Our Main office is both Nikon and Olympus Authorized for Service and Warranty Repairs.

Bi-Annual Service is recommended to reduce downtime and keep your equipment in exceptional shape. We often service brands such as Nikon, Olympus, Zeiss, Leica, Motic, Accu-Scope Unitron, LW Scientific, Swift, and National. Most replacement parts are stocked to reduce downtime.

Many microscope manufacturers quit producing replacement parts after a set period of time. Most of which are still great microscopes that just need attention to keep running. Objectives can be rebuilt or replaced, Illumination systems can be upgraded to LED or restored to original. If a gear is damaged, and a replacement is not available, we can cast a replacement. Optical alignment can be done to reduce double image or image shift as magnification increases.

New microscopes are available from Accu-Scope, Unitron, The Motic Group and Meiji. We also offer Jenoptik, Lumenera, Excelis and Moticam Cameras. Are you interested in a Microscope or Camera brand that we do not carry? Other brands are available through our dealer network. Every microscope user has different demands, lets start a conversation about your needs and how we can get there.



Main Office

Faribault, Minnesota

Southern Office

Keller, Texas

Satelite Office

Huntington Beach, CA

Kansas City, KS

Dayton, OH

Harrisburg, PA

Tampa, FL




Pay Online
Credit Card
Check
PO

BMIC

BLAIS MICROSCOPE COMPANY

Microscopy Solutions

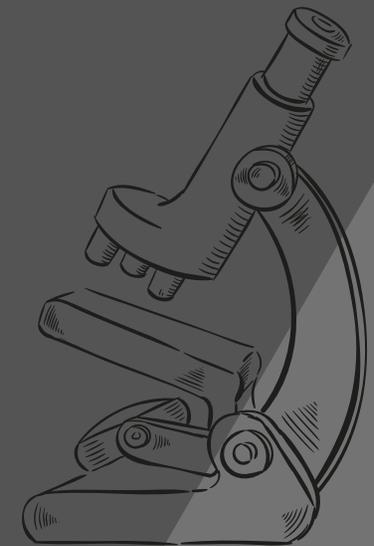

Sales Service Repair

(800) 223-2259

info@blaismicroscope.com

www.blaismicroscope.com

Microscope Service

Objectives, Condenser, Slide Holder and Accessories removed




Microscope and Components scrubbed and cleaned



Objectives are cleaned and reinstalled. Shimed to parfocal as needed



Eyepieces and View Windows in Head cleaned

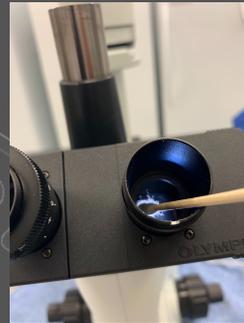

Slide Holder is reinstalled and adjusted to ensure proper scanning

Exposed Components lubricated as need

Condenser is reinstalled and centered (Koehler illumination when applicable)

Electrical System checked for proper functionality



Check for proper bulb



Microscope is tested with a slide, any issues will then be addressed



Get your Microscope Serviced today! Price includes travel when scheduling on route with our other customers.